

Cycle Control

By K & M Products

K & M Products Cycle Control automates Chatillon® LTCM-100 & LTCM-500 Series Testers



Cycle Control conveniently plugs into the tester footswitch jack, enhancing the economical LTCM Series tester by providing 3 additional modes:

- In **Manual Mode** the tester will travel in the desired direction and stop at the preset mechanical limit.
- In **Single Cycle Mode** the tester will travel in the desired beginning direction and upon reaching the preset limit will reverse direction, travel to the preset opposing limit and stop.
- In **Cycle Mode** the tester will cycle continuously between preset upper/lower limits.
- **Cycle Control** also functions with the load limit feature on Chatillon® DF Series Digital Force Gages.

Specifications:

- Requires Chatillon® LTCM-100 or LTCM-500 Tester with footswitch jack.
- 110V 50/60 Hz operation.
- Operating temperature range: 40°F to 100°F (5°C to 45°C).
- Dimension 4-3/8" x 4-3/8" x 4".

Cycle Control installation is simple and requires no modifications to the tester.

Cycle Control works with any Chatillon® LTCM-100 or LTCM-500 Tester having footswitch jack capability.

Features:


- Manual Control
- Single Cycle
- Cycle
- Lighted up/down travel indicator
- Fast & Easy Installation
- 1 Year Limited Warranty

Includes:

- Power Cable
- Footswitch Jack Cable
- User Manual

Press and Release operation

- **Cycle Control** eliminates the need to hold the tester up/down switch while the test is running!
- **Cycle Control** has integral lamps in the up/down switch indicating direction of travel.
- **Cycle Control** automates your test whether using Load Limit or Mechanical Limit Tests.
- **Cycle Control** automatically selects return direction in single cycle and cycle modes.

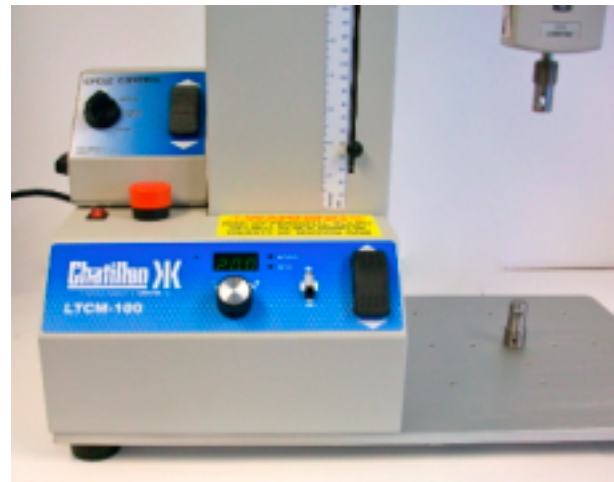

K & M Products

A department of S. A. Meier Co. of Milwaukee, Inc.
230 James Street, A-1 / Wales, WI 53183 / 262 968-4950 / sales @samco-inc.com

03/08